

CDN Resource Laboratories Ltd.

#2, 20148 – 102nd Avenue, Langley, B.C., Canada, V1M 4B4, 604-882-8422, Fax: 604-882-8466 (www.cdnlabs.com)

STANDARD REFERENCE MATERIAL: CDN-GS-5Q

Recommended values and the “Between Lab” Two Standard Deviations

<i>Gold</i>	<i>5.59 g/t ± 0.35 g/t</i>	<i>30g FA, instrumental</i>	<i>Certified value</i>
<i>Gold</i>	<i>5.64 g/t ± 0.35 g/t</i>	<i>30g FA, gravimetric</i>	<i>Certified value</i>
<i>Silver</i>	<i>60.3 g/t ± 3.9 g/t</i>	<i>4-acid, instrumental</i>	<i>Certified value</i>

PREPARED BY: CDN Resource Laboratories Ltd.

CERTIFIED BY: Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia

INDEPENDENT GEOCHEMIST: Dr. Barry Smee., Ph.D., P. Geo.

DATE OF CERTIFICATION: February 15, 2014

ORIGIN OF REFERENCE MATERIAL:

Standard CDN-GS-5Q was prepared using 740 kg of granitic rock blended with 60 kg of a high grade Au-Ag ore.

METHOD OF PREPARATION:

Reject ore material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 material was mixed for 5 days in a double-cone blender. Splits were taken and sent to 15 commercial laboratories for round robin assaying.

Approximate chemical composition (by whole rock analysis) is as follows:

	Percent		Percent
SiO ₂	58.9	MgO	2.6
Al ₂ O ₃	15.9	K ₂ O	1.5
Fe ₂ O ₃	8.4	TiO ₂	0.5
CaO	5.8	LOI	1.5
Na ₂ O	3.1	S	0.3

Statistical Procedures:

The final limits were calculated after first determining if all data was compatible within a spread normally expected for similar analytical methods done by reputable laboratories. Data from any one laboratory was removed from further calculations when the mean of all analyses from that laboratory failed a t test of the global means of the other laboratories. The means and standard deviations were calculated using all remaining data. Any analysis that fell outside of the mean ± 2 standard deviations was removed from the ensuing data base. The mean and standard deviations were again calculated using the remaining data. This method is different from that used by Government agencies in that the actual “between-laboratory” standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.

Results from round-robin assaying are displayed on the following page.

STANDARD REFERENCE MATERIAL CDN-GS-5Q

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
Instrumental	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
GS-5Q-1	5.35	5.71	5.62	5.38	5.48	5.74	5.50	5.71	5.09	5.54	5.23	5.62	5.89		5.64
GS-5Q-2	5.62	5.75	5.78	5.84	5.81	5.60	5.50	5.72	5.35	5.80	5.39	5.70	5.30		5.61
GS-5Q-3	5.63	5.58	5.86	5.32	5.64	5.72	5.40	5.61	5.63	5.82	5.44	5.43	5.70		5.67
GS-5Q-4	5.61	5.80	5.82	5.77	5.48	5.72	6.18	6.09	5.24	5.38	5.14	5.37	5.79		5.71
GS-5Q-5	5.74	5.80	5.74	5.72	5.58	5.69	6.21	5.78	5.37	5.35	5.21	5.24	5.47		5.58
GS-5Q-6	5.69	5.61	5.66	5.69	5.60	5.70	6.83	5.76	5.83	5.58	5.41	5.60	5.48		5.66
GS-5Q-7	5.64	5.70	5.68	5.79	5.59	5.73	5.41	5.88	5.41	5.75	5.29	4.72	5.31		5.73
GS-5Q-8	5.75	5.68	5.50	5.61	5.72	5.64	5.32	5.61	5.33	5.17	5.35	5.85	5.51		5.72
GS-5Q-9	5.62	5.60	5.90	5.57	5.71	5.64	5.28	5.56	5.39	5.50	5.29	6.14	5.58		5.61
GS-5Q-10	5.90	5.70	5.69	5.65	5.43	5.72	5.20	5.76	5.53	5.13	5.53	5.29	5.21		5.70
Mean	5.66	5.69	5.72	5.63	5.60	5.69	5.68	5.75	5.42	5.50	5.33	5.50	5.52		5.66
Std. Dev'n	0.1398	0.0765	0.1197	0.1695	0.1197	0.0471	0.5364	0.1532	0.2062	0.2462	0.1184	0.3856	0.2207		0.0525
%RSD	2.47	1.34	2.09	3.01	2.14	0.83	9.44	2.67	3.81	4.47	2.22	7.02	4.00		0.93
Gravimetric	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
GS-5Q-1	5.51	5.80	5.72	5.82	5.78	5.99	5.47	5.73	5.35	5.72	5.61	5.33	5.77	5.85	5.43
GS-5Q-2	5.66	5.70	5.62	5.66	5.74	5.58	5.30	5.93	5.45	5.49	5.79	5.57	5.44	5.57	5.48
GS-5Q-3	5.48	5.60	5.73	5.35	5.62	5.99	5.67	5.94	5.59	5.76	5.69	5.57	5.61	5.71	5.49
GS-5Q-4	5.76	5.60	5.83	5.50	5.63	5.68	5.83	5.89	5.66	5.75	6.16	5.46	5.59	5.89	5.53
GS-5Q-5	5.76	5.70	5.49	5.58	5.48	5.58	5.80	5.79	5.49	5.68	5.61	5.79	5.48	5.63	5.42
GS-5Q-6	5.74	5.70	5.62	5.86	5.48	6.09	6.13	6.14	5.38	5.55	5.64	5.51	5.57	5.53	5.41
GS-5Q-7	5.59	5.70	5.36	5.57	5.78	6.09	5.23	6.09	5.14	5.51	5.95	5.44	5.72	5.80	5.46
GS-5Q-8	5.69	5.90	5.54	5.84	5.49	5.34	5.87	5.98	5.55	5.51	6.04	5.69	5.46	5.80	5.56
GS-5Q-9	5.70	5.40	5.59	5.62	5.59	5.83	5.67	6.08	5.59	5.55	5.82	5.46	5.70	5.88	5.59
GS-5Q-10	5.70	5.60	5.61	5.76	5.59	6.10	5.73	6.05	5.59	5.68	5.72	5.49	5.73	5.51	5.47
Mean	5.66	5.67	5.61	5.66	5.62	5.83	5.67	5.96	5.48	5.62	5.80	5.53	5.61	5.72	5.48
Std. Dev'n	0.100	0.134	0.132	0.164	0.117	0.268	0.272	0.133	0.155	0.108	0.190	0.132	0.120	0.147	0.060
%RSD	1.77	2.36	2.36	2.91	2.09	4.60	4.81	2.24	2.83	1.92	3.28	2.38	2.15	2.57	1.09
	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t
GS-5Q-1	59.6	60.7	61	59.2	63	60.5	62.7	57.1	58.4	62	52.6	57	63.1	56.0	58.4
GS-5Q-2	61.4	58.5	60	58.8	63	60.1	62.5	58.0	57.6	64	53.9	60	61.3	57.3	58.3
GS-5Q-3	59.7	59.6	62	59.5	63	62.0	62.6	57.4	58.5	63	52.8	61	62.6	58.8	59.9
GS-5Q-4	61.0	59.0	60	58.7	65	60.8	62.1	58.5	62.4	61	53.6	59	61.3	57.6	57.6
GS-5Q-5	58.4	60.5	61	58.7	64	61.1	62.1	58.5	57.3	62	53.4	62	59.9	57.4	58.1
GS-5Q-6	60.7	60.0	62	58.1	64	63.1	62.4	58.7	60.3	61	54.4	61	59.5	58.4	57.5
GS-5Q-7	61.6	60.0	60	59.2	64	60.6	61.9	58.6	57.2	63	52.3	60	61.4	57.3	60.3
GS-5Q-8	60.1	58.2	62	58.8	63	60.6	60.8	59.2	57.3	64	55.0	60	62.1	57.1	57.9
GS-5Q-9	60.8	60.1	61	58.0	63	60.2	63.3	58.6	63.4	62	55.7	61	63.4	57.1	59.5
GS-5Q-10	59.2	61.0	60	58.4	65	58.4	63.8	58.1	60.4	65	53.8	59	64.1	58.1	57.4
Mean	60.3	59.8	60.9	58.7	63.7	60.7	62.4	58.3	59.3	62.7	53.8	60.0	61.9	57.5	58.5
Std. Dev'n	1.027	0.932	0.876	0.481	0.823	1.228	0.809	0.633	2.245	1.337	1.069	1.414	1.483	0.784	1.043
%RSD	1.70	1.56	1.44	0.82	1.29	2.02	1.30	1.09	3.79	2.13	1.99	2.36	2.40	1.36	1.78

Notes:

- 1) **Ag data from laboratory 11 was excluded from the calculations for failing the t test.**
- 2) **Lab 14 was unable to supply Au data with an instrumental finish.**

STANDARD REFERENCE MATERIAL CDN-GS-5Q

Participating Laboratories:


(not in same order as listed in table of results)

Acme Analytical Laboratories Ltd., Vancouver, B.C., Canada
Acme Analytical Laboratories Ltd, Santiago, Chile
Activation Laboratories, Ancaster, Ontario, Canada
Activation Laboratories, Kamloops, B.C., Canada
Activation Laboratories, Thunder Bay, Ontario, Canada
ALS Canada, North Vancouver, B.C., Canada
ALS, Loughrea, Ireland
American Assay Lab., Nevada, USA
Alex Stewart Assayers, Mendoza, Argentina
Certimin, Lima, Peru
Met-Solve Analytical Services, Langley, B.C., Canada
Skyline Assayers and Laboratories, Arizona, USA
SGS – Lima, Peru
SGS, Vancouver, B.C., Canada
TSL Laboratories Ltd., Saskatoon, SK, Canada

Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. or Barry Smee accept no liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by

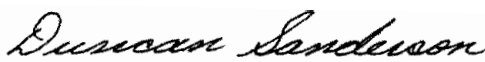

Duncan Sanderson, Certified Assayer of B.C.

Geochemist



Dr. Barry Smee, Ph.D., P. Geo.