

CDN Resource Laboratories Ltd.

#2, 20148 - 102nd Avenue, Langley, B.C., Canada, V1M 4B4, Ph: 604-882-8422 Fax: 604-882-8466
(www.cdnlabs.com)

STANDARD REFERENCE MATERIAL: CDN-BL-7

Recommended values:

Gold concentration: < 0.01 g/t

Platinum concentration: < 0.01 g/t

Palladium concentration: < 0.01 g/t

PREPARED BY: CDN Resource Laboratories Ltd.
CERTIFIED BY: Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia
INDEPENDENT GEOCHEMIST: Dr. Barry Smee., Ph. D., P. Geo.
DATE OF CERTIFICATION: June 22, 2010

ORIGIN OF REFERENCE MATERIAL:

Standard CDN-BL-7 was prepared using a blank granitic material.

METHOD OF PREPARATION:

The granitic material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 (<53 micron) material was mixed for 5 days in a double-cone blender. Splits were taken and sent to 6 commercial laboratories for round robin assaying. Round robin results are displayed on the next page.

APPROXIMATE CHEMICAL COMPOSITION (whole rock analysis):

| | Percent | | | Percent |
|--------------------------------|---------|--|-------------------|---------|
| SiO ₂ | 73.2 | | Na ₂ O | 2.9 |
| Al ₂ O ₃ | 10.7 | | MgO | 1.6 |
| Fe ₂ O ₃ | 4.0 | | K ₂ O | 1.3 |
| CaO | 2.6 | | TiO ₂ | 0.4 |
| MnO | 0.1 | | LOI | 1.5 |
| C | 0.1 | | S | 0.1 |

Statistical Procedures: There was no statistical analysis performed on the data.

Participating Laboratories: (not in same order as table of assays)

Acme Analytical Laboratories Ltd., Vancouver
Assayers Canada Ltd., Vancouver
ALS Chemex Laboratories, North Vancouver
Eco Tech Laboratory Ltd., Kamloops, B.C., Canada
International Plasma Labs, Richmond, B.C., Canada
TSL Laboratories, Saskatoon

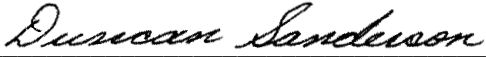
Assay Procedure: assays were fire assay, AA or ICP finish on 30g samples.
Assay results on the following page

| | Lab 1 | Lab 2 | Lab 3 | Lab 4 | Lab 5 | Lab 6 |
|-------------|--------|--------|--------|--------|--------|--------|
| Sample | Au ppb | Au ppb | Au ppb | Au ppb | Au ppb | Au ppb |
| CDN-BL-7-1 | 6 | <1 | <5 | <5 | 8 | <5 |
| CDN-BL-7-2 | 3 | 1 | <5 | <5 | <5 | <5 |
| CDN-BL-7-3 | <2 | 5 | <5 | <5 | <5 | <5 |
| CDN-BL-7-4 | <2 | 2 | <5 | <5 | 8 | <5 |
| CDN-BL-7-5 | 5 | 2 | <5 | <5 | 5 | <5 |
| CDN-BL-7-6 | <2 | 1 | <5 | <5 | <5 | <5 |
| CDN-BL-7-7 | <2 | 2 | <5 | <5 | 11 | <5 |
| CDN-BL-7-8 | 3 | 1 | <5 | <5 | 5 | <5 |
| CDN-BL-7-9 | <2 | 1 | <5 | <5 | 8 | <5 |
| CDN-BL-7-10 | <2 | 2 | <5 | <5 | <5 | <5 |
| | | | | | | |
| | Pt ppb | Pt ppb | Pt ppb | Pt ppb | Pt ppb | Pt ppb |
| CDN-BL-7-1 | <3 | 5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-2 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-3 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-4 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-5 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-6 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-7 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-8 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-9 | <3 | <5 | <5 | <5 | <5 | <10 |
| CDN-BL-7-10 | 4 | <5 | <5 | <5 | <5 | <10 |
| | | | | | | |
| | Pd ppb | Pd ppb | Pd ppb | Pd ppb | Pd ppb | Pd ppb |
| CDN-BL-7-1 | <2 | 1 | <5 | <5 | 10 | <5 |
| CDN-BL-7-2 | <2 | 2 | <5 | <5 | <5 | <5 |
| CDN-BL-7-3 | <2 | 2 | <5 | <5 | 14 | <5 |
| CDN-BL-7-4 | <2 | 1 | <5 | <5 | 7 | <5 |
| CDN-BL-7-5 | <2 | 1 | <5 | <5 | <5 | <5 |
| CDN-BL-7-6 | <2 | 1 | <5 | <5 | 8 | <5 |
| CDN-BL-7-7 | <2 | 2 | <5 | <5 | <5 | <5 |
| CDN-BL-7-8 | <2 | 1 | <5 | <5 | 6 | <5 |
| CDN-BL-7-9 | <2 | 1 | <5 | <5 | 9 | <5 |
| CDN-BL-7-10 | <2 | 1 | <5 | <5 | 8 | <5 |


Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. nor Barry Smee accept any liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by


 Duncan Sanderson, Certified Assayer of B.C.

Geochemist


 Dr. Barry Smee, Ph.D., P. Geo.