

CDN Resource Laboratories Ltd.

#2, 20148 – 102nd Avenue, Langley, B.C., Canada, V1M 4B4, 604-882-8422, Fax: 604-882-8466 (www.cdnlabs.com)

ORE REFERENCE STANDARD: CDN-CM-8

Recommended values and the “Between Lab” Two Standard Deviations

Gold: *0.91 ± 0.11 g/t*
Copper: *0.364 ± 0.024 %*
Molybdenum: *0.0160 ± 0.0014 %*

PREPARED BY: CDN Resource Laboratories Ltd.
CERTIFIED BY: Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia
INDEPENDENT GEOCHEMIST: Dr. Barry Smee., Ph.D., P. Geo.
DATE OF CERTIFICATION: June 14, 2010

ORIGIN OF REFERENCE MATERIAL:

Standard CDN-CM-8 was prepared using a North American calc-alkalic copper-gold-molybdenum porphyry ore. It is derived from altered granodiorite, mafic to intermediate volcanic and volcanoclastic sedimentary rocks. Mineralization is principally pyrite, chalcopyrite and molybdenite that occurs in veins, stockworks and disseminations. 792 kg of this ore was blended with 2 kg of a high grade gold ore.

METHOD OF PREPARATION:

Reject ore material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 material was mixed for 5 days in a double-cone blender. Splits were taken and sent to 15 commercial laboratories for round robin assaying.

Approximate chemical composition (by whole rock analysis) is as follows:

	Percent			Percent
SiO ₂	60.7		MgO	1.5
Al ₂ O ₃	15.0		K ₂ O	6.9
Fe ₂ O ₃	7.5		TiO ₂	0.4
CaO	0.7		LOI	4.5
Na ₂ O	1.5		S	2.6
C	0.3			

Statistical Procedures:

The final limits were calculated after first determining if all data was compatible within a spread normally expected for similar analytical methods done by reputable laboratories. Data from any one laboratory was removed from further calculations when the mean of all analyses from that laboratory failed a t test of the global means of the other laboratories. The means and standard deviations were calculated using all remaining data. Any analysis that fell outside of the mean ±2 standard deviations was removed from the ensuing data base. The mean and standard deviations were again calculated using the remaining data. This method is different from that used by Government agencies in that the actual “between-laboratory” standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.

Results from round-robin assaying are displayed on the following page.

STANDARD REFERENCE MATERIAL CDN-CM-8

Assay Procedures: **Au:** Fire assay pre-concentration, AA or ICP finish (30g sub-sample).

Cu, Mo: 4-acid digestion, AA or ICP finish.

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
SAMPLE	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
CDN-CM-8-1	0.85	0.810	0.97	0.94	0.81	0.770	1.010	0.985	0.95	0.932	0.96	0.93	0.909	0.947	0.90
CDN-CM-8-2	0.82	0.896	0.96	0.89	0.90	0.835	0.962	0.914	0.91	0.842	0.96	0.86	0.940	0.844	0.84
CDN-CM-8-3	0.94	0.925	1.00	0.88	0.82	0.875	0.972	0.991	0.92	0.802	0.92	0.92	0.908	0.920	0.91
CDN-CM-8-4	0.89	0.983	0.96	0.86	0.79	0.877	0.961	0.979	0.94	0.850	0.99	0.95	0.910	1.031	0.92
CDN-CM-8-5	0.88	0.865	0.93	0.88	0.82	0.830	0.952	0.920	0.92	0.765	0.89	0.90	0.903	0.851	0.83
CDN-CM-8-6	0.82	0.827	0.92	0.90	0.97	0.822	0.990	0.955	0.96	0.807	0.89	0.82	0.888	0.870	0.92
CDN-CM-8-7	0.94	0.886	0.96	0.81	0.80	0.869	1.090	1.000	0.88	0.939	0.86	1.01	0.906	0.810	0.91
CDN-CM-8-8	0.85	0.993	0.97	0.78	0.84	0.847	0.989	0.937	0.96	0.821	0.92	1.00	0.965	0.879	0.89
CDN-CM-8-9	0.88	0.949	0.95	0.86	0.91	0.754	1.010	0.915	0.95	0.822	0.96	0.93	0.897	0.867	0.97
CDN-CM-8-10	0.89	1.025	1.00	0.88	0.89	0.839	0.984	0.922	0.93	0.828	0.92	0.86	0.919	0.975	0.88
Mean	0.876	0.916	0.962	0.868	0.855	0.832	0.992	0.952	0.932	0.841	0.927	0.918	0.915	0.899	0.897
Std. Dev'n	0.0425	0.0719	0.0257	0.0452	0.0591	0.0416	0.0397	0.0343	0.0253	0.0551	0.0403	0.0607	0.0224	0.068	0.04
%RSD	4.85	7.85	2.68	5.20	6.91	5.00	4.00	3.61	2.71	6.56	4.35	6.61	2.45	7.55	4.52
	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %	Cu %
CDN-CM-8-1	0.358	0.360	0.361	0.37		0.374	0.317	0.333	0.355		0.37	0.357	0.373	0.367	0.364
CDN-CM-8-2	0.362	0.357	0.362	0.37		0.373	0.340	0.327	0.344		0.37	0.353	0.368	0.360	0.377
CDN-CM-8-3	0.354	0.363	0.363	0.36		0.375	0.330	0.326	0.346		0.37	0.360	0.369	0.363	0.380
CDN-CM-8-4	0.361	0.368	0.360	0.38		0.370	0.328	0.327	0.355		0.37	0.356	0.370	0.372	0.365
CDN-CM-8-5	0.362	0.375	0.366	0.37		0.364	0.324	0.326	0.354		0.37	0.363	0.371	0.373	0.370
CDN-CM-8-6	0.362	0.356	0.361	0.38		0.373	0.323	0.329	0.344		0.37	0.353	0.376	0.365	0.370
CDN-CM-8-7	0.360	0.347	0.361	0.37		0.374	0.319	0.334	0.355		0.38	0.354	0.377	0.374	0.368
CDN-CM-8-8	0.369	0.365	0.368	0.38		0.372	0.340	0.333	0.350		0.37	0.362	0.378	0.360	0.382
CDN-CM-8-9	0.374	0.370	0.368	0.37		0.368	0.350	0.341	0.335		0.38	0.362	0.372	0.368	0.368
CDN-CM-8-10	0.374	0.372	0.369	0.38		0.368	0.329	0.341	0.341		0.38	0.358	0.383	0.368	0.377
Mean	0.364	0.363	0.364	0.373		0.371	0.330	0.332	0.348		0.373	0.358	0.374	0.367	0.372
Std. Dev'n	0.0066	0.0085	0.0035	0.0067		0.0035	0.0104	0.0058	0.0070		0.0048	0.0038	0.0047	0.0051	0.0064
%RSD	1.83	2.34	0.96	1.81		0.95	3.16	1.74	2.01		1.30	1.07	1.26	1.38	1.72
	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %	Mo %
CDN-CM-8-1	0.016	0.0163	0.016	0.017		0.017	0.012	0.014	0.0144		0.016	0.0161	0.0151	0.0163	0.0160
CDN-CM-8-2	0.016	0.0153	0.015	0.017		0.016	0.013	0.014	0.0131		0.016	0.0155	0.0162	0.0157	0.0165
CDN-CM-8-3	0.015	0.0151	0.016	0.017		0.017	0.013	0.015	0.0138		0.016	0.0167	0.0159	0.0164	0.0165
CDN-CM-8-4	0.016	0.0159	0.015	0.016		0.017	0.013	0.015	0.0152		0.017	0.0161	0.0157	0.0159	0.0175
CDN-CM-8-5	0.016	0.0160	0.016	0.016		0.017	0.013	0.015	0.0143		0.017	0.0163	0.0152	0.0158	0.0165
CDN-CM-8-6	0.015	0.0149	0.016	0.016		0.017	0.014	0.014	0.0146		0.017	0.0161	0.0162	0.0159	0.0165
CDN-CM-8-7	0.015	0.0145	0.016	0.016		0.017	0.012	0.014	0.0157		0.017	0.0171	0.0163	0.0165	0.0170
CDN-CM-8-8	0.016	0.0156	0.015	0.016		0.017	0.013	0.015	0.0153		0.016	0.0164	0.0152	0.0161	0.0170
CDN-CM-8-9	0.016	0.0155	0.015	0.016		0.017	0.013	0.015	0.0141		0.017	0.0163	0.0153	0.0165	0.0170
CDN-CM-8-10	0.016	0.0156	0.016	0.016		0.017	0.014	0.013	0.0143		0.016	0.0158	0.0154	0.0160	0.0170
Mean	0.0157	0.0155	0.0156	0.0163		0.0169	0.0130	0.0144	0.0145		0.0165	0.0162	0.0157	0.0161	0.0168
Std. Dev'n	0.0005	0.0005	0.0005	0.0005		0.0003	0.0007	0.0007	0.0008		0.0005	0.0004	0.0005	0.0003	0.0004
%RSD	3.08	3.49	3.31	2.96		1.87	5.13	4.86	5.29		3.19	2.74	3.00	1.84	2.54

Note: "Mo" data from laboratory 7 was excluded from the calculations for failing the t test.
Laboratories 5 and 10 did not report Cu nor Mo data.

STANDARD REFERENCE MATERIAL CDN-CM-8

Participating Laboratories:


(not in same order as listed in table of results)

Acme Analytical Laboratories Ltd., Vancouver, B.C., Canada
Activation Laboratories, Ancaster, Ontario, Canada
Activation Laboratories, Thunder Bay, Ontario, Canada
ALS Chemex, North Vancouver, B.C., Canada
ALS Chemex, Nevada, USA
Assayers Canada Ltd., Vancouver, B.C., Canada
Bourlamaque Assay Laboratories, Quebec, Canada
Eco Tech Laboratory Ltd., Kamloops, B.C., Canada
Genalysis, Perth, Australia
Labtium Inc., Finland
Omac, Ireland
SGS, Toronto, Canada
Skyline, Tucson, Arizona
TSL Laboratories Ltd., Saskatoon, SK, Canada
Ultra Trace, Perth, Australia

Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. or Barry Smee accept no liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by

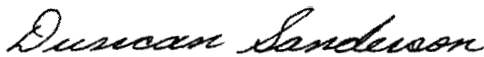

Duncan Sanderson, Certified Assayer of B.C.

Geochemist



Dr. Barry Smee, Ph.D., P. Geo.