

CDN Resource Laboratories Ltd.

10945-B River Road, Delta, B.C., Canada, V4C 2R8, 604-540-2233, Fax: 604-588-3960 (www.cdnlabs.com)

ORE REFERENCE STANDARD: CDN-FCM-3

Recommended values and the "Between Lab" Two Standard Deviations

Gold 0.40 ± 0.07 g/t
Silver 23.6 ± 3.3 g/t
Copper 0.291 ± 0.020 %
Lead 0.152 ± 0.014 %
Zinc 0.543 ± 0.032 %

PREPARED BY: CDN Resource Laboratories Ltd.
CERTIFIED BY: Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia
INDEPENDENT GEOCHEMIST: Dr. Barry Smee., Ph.D., P. Geo.
DATE OF CERTIFICATION: February 20, 2006

METHOD OF PREPARATION:

Reject ore material was dried, crushed, pulverized and then passed through a 200 mesh screen. The +200 material was discarded. The -200 material was mixed for 5 days in a rotary mixer. Splits were taken and sent to twelve laboratories for round robin assaying. The material has been packaged in nominal 100g lots in tin-top kraft bags which have been individually vacuum-sealed in polyethylene bags.

ORIGIN OF REFERENCE MATERIAL:

The ore was supplied by Hunter Dickinson (Farallon) from their Campo Morado property in Mexico. The Campo Morado precious-metal-bearing, volcanogenic massive sulphide deposits occur in a lower Cretaceous bimodal, calc-alkaline volcanic sequence. Most deposits occur in the upper part of a sequence of felsic flows and heterolithic volcanoclastic rocks or at its contact with overlying chert and argillite. Gold, silver, zinc, and lead are associated with pyrite, quartz, ankerite, sphalerite, chalcopyrite and galena, with minor tennantite-freibergite, arsenopyrite, and pyrrotite. Standard CDN-FCM-3 was made by combining 145 kg of Farallon material with 218 kg of blank granitic ore.

Approximate chemical composition is as follows:

	Percent		Percent
SiO ₂	48.2	MgO	1.0
Al ₂ O ₃	8.6	K ₂ O	1.1
Fe ₂ O ₃	22.3	TiO ₂	0.2
CaO	2.2	LOI	11.5
Na ₂ O	2.7		

Statistical Procedures:

The mean and standard deviation for all data was calculated. Outliers were defined as samples beyond the mean ± 2 Standard Deviations from all data. These outliers were removed from the data and a new mean and standard deviation was determined. The Au data from one laboratory and the Ag data from another laboratory were excluded as they did not pass the "t" test. This method is different from that used by Government agencies in that the actual "between-laboratory" standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.

Results from round-robin assaying are presented on subsequent pages:

Assay Procedures:

Au: Fire assay pre-concentration, AA or ICP finish (10g sub-sample).
Ag, Cu, Pb, Zn: 4-acid digestion, AA or ICP finish.

STANDARD REFERENCE MATERIAL CDN-FCM-3

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12
	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
FCM-3-1	0.38	0.44	0.47	0.41	0.48	0.39	0.42	0.43	0.44	0.40	0.35	0.35
FCM-3-2	0.39	0.42	0.47	0.39	0.44	0.41	0.45	0.42	0.43	0.40	0.36	0.34
FCM-3-3	0.38	0.40	0.53	0.42	0.43	0.41	0.46	0.41	0.42	0.37	0.34	0.36
FCM-3-4	0.45	0.39	0.42	0.37	0.43	0.38	0.46	0.48	0.42	0.38	0.34	0.33
FCM-3-5	0.43	0.42	0.47	0.38	0.42	0.42	0.43	0.39	0.43	0.38	0.36	0.37
FCM-3-6	0.43	0.41	0.40	0.38	0.45	0.38	0.43	0.40	0.42	0.37	0.34	0.37
FCM-3-7	0.39	0.42	0.42	0.39	0.43	0.38	0.44	0.38	0.42	0.40	0.37	0.36
FCM-3-8	0.37	0.42	0.49	0.38	0.44	0.40	0.45	0.38	0.42	0.38	0.35	0.36
FCM-3-9	0.38	0.39	0.47	0.41	0.41	0.41	0.45	0.38	0.44	0.46	0.36	0.35
FCM-3-10	0.44	0.39	0.41	0.38	0.42	0.39	0.45	0.42	0.42	0.40	0.35	0.37
Mean	0.40	0.41	0.46	0.39	0.44	0.40	0.44	0.41	0.43	0.39	0.35	0.36
Std. Devn.	0.0299	0.0170	0.0412	0.0169	0.0196	0.0149	0.0127	0.0311	0.0084	0.0263	0.0091	0.0135
% RSD	7.40	4.15	9.05	4.36	4.50	3.76	2.87	7.60	1.98	6.68	2.59	3.79
	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t
FCM-3-1	21.9	20.9	27	24	25	21	25.7	23.5	23.9	23	23.1	22.6
FCM-3-2	22.3	20.8	26	25	26	22	25.7	23.7	23.9	25	23.2	23.2
FCM-3-3	21.6	20.7	30	25	26	23	26.2	24.3	23.7	24	22.5	22.9
FCM-3-4	21.8	20.9	27	24	26	22	25.9	22.8	24.4	23	22.6	22.6
FCM-3-5	21.6	21.0	28	24	26	23	26.2		23.5	23	21.3	24.0
FCM-3-6	22.1	20.9	29	25	25	23	25.1	24.7	23.8	25	21.2	22.7
FCM-3-7	21.9	20.7	25	23	26	23	25.4	24.2	23.8	23	22.7	22.7
FCM-3-8	22.3	20.8	26	25	26	23	25.9	23.9	23.7	23	22.3	22.6
FCM-3-9	21.7	20.9	27	23	25	22	25.4	25.7	23.8	23	22.4	22.7
FCM-3-10	21.6	20.9	28	24	26	23	26.3	24.5	23.8	23	22.6	22.9
Mean	21.9	20.9	27.3	24.2	25.7	22.5	25.8	24.1	23.8	23.6	22.4	22.9
Std. Devn.	0.2741	0.0972	1.4944	0.7888	0.4830	0.7071	0.3967	0.8187	0.2312	0.6835	0.6641	0.4332
% RSD	1.25	0.47	5.47	3.26	1.88	3.14	1.54	3.39	0.97	2.90	2.97	1.89

STANDARD REFERENCE MATERIAL CDN-FCM-3

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12
	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu
FCM-3-1	0.278	0.28	0.301	0.295	0.280	0.299	0.299	0.298	0.28	0.291	0.30	0.290
FCM-3-2	0.274	0.28	0.301	0.282	0.283	0.304	0.302	0.301	0.28	0.292	0.30	0.293
FCM-3-3	0.280	0.28	0.310	0.283	0.282	0.306	0.291	0.300	0.28	0.292	0.30	0.292
FCM-3-4	0.272	0.29	0.304	0.291	0.283	0.309	0.297	0.296	0.28	0.294	0.31	0.291
FCM-3-5	0.275	0.28	0.302	0.283	0.283	0.312	0.304		0.28	0.290	0.31	0.293
FCM-3-6	0.276	0.28	0.307	0.276	0.280	0.306	0.292	0.299	0.29	0.290	0.31	0.290
FCM-3-7	0.279	0.29	0.302	0.280	0.281	0.300	0.290	0.295	0.29	0.294	0.30	0.292
FCM-3-8	0.278	0.28	0.300	0.283	0.284	0.299	0.294	0.300	0.28	0.294	0.31	0.291
FCM-3-9	0.274	0.28	0.304	0.287	0.282	0.303	0.291	0.296	0.28	0.289	0.30	0.290
FCM-3-10	0.276	0.28	0.299	0.286	0.281	0.305	0.293	0.300	0.28	0.291	0.31	0.291
Mean	0.276	0.282	0.303	0.285	0.282	0.304	0.295	0.298	0.282	0.292	0.305	0.291
Std. Devn.	0.0025	0.0042	0.0034	0.0054	0.0014	0.0043	0.0049	0.0022	0.0042	0.0018	0.0053	0.0012
% RSD	0.92	1.50	1.11	1.91	0.49	1.40	1.67	0.73	1.50	0.62	1.73	0.40
	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb
FCM-3-1	0.16	0.14	0.16	0.158	0.151	0.146	0.164	0.156	0.16	0.150	0.15	0.147
FCM-3-2	0.16	0.13	0.17	0.147	0.151	0.147	0.164	0.156	0.15	0.150	0.15	0.146
FCM-3-3	0.16	0.14	0.17	0.149	0.151	0.149	0.168	0.155	0.15	0.153	0.15	0.146
FCM-3-4	0.16	0.14	0.16	0.152	0.151	0.151	0.166	0.155	0.15	0.154	0.16	0.146
FCM-3-5	0.16	0.14	0.16	0.150	0.152	0.154	0.169		0.15	0.150	0.16	0.150
FCM-3-6	0.16	0.14	0.17	0.144	0.150	0.148	0.166	0.154	0.15	0.149	0.16	0.145
FCM-3-7	0.17	0.14	0.17	0.148	0.150	0.146	0.166	0.157	0.15	0.150	0.15	0.148
FCM-3-8	0.15	0.14	0.17	0.147	0.151	0.144	0.165	0.157	0.15	0.153	0.15	0.146
FCM-3-9	0.15	0.14	0.16	0.151	0.151	0.148	0.166	0.154	0.15	0.150	0.15	0.146
FCM-3-10	0.15	0.14	0.16	0.150	0.150	0.150	0.163	0.155	0.15	0.150	0.15	0.145
Mean	0.158	0.139	0.165	0.150	0.151	0.148	0.166	0.155	0.151	0.151	0.153	0.147
Std. Devn.	0.0063	0.0032	0.0053	0.0037	0.0006	0.0029	0.0018	0.0011	0.0032	0.0018	0.0048	0.0015
% RSD	4.00	2.28	3.19	2.51	0.42	1.93	1.10	0.73	2.09	1.17	3.16	1.03
	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn
FCM-3-1	0.52	0.50	0.56	0.553	0.530	0.553	0.573	0.555	0.54	0.546	0.53	0.523
FCM-3-2	0.53	0.51	0.55	0.544	0.534	0.538	0.555	0.556	0.55	0.555	0.54	0.528
FCM-3-3	0.53	0.52	0.57	0.541	0.530	0.550	0.560	0.554	0.55	0.563	0.55	0.526
FCM-3-4	0.53	0.52	0.56	0.559	0.533	0.558	0.566	0.561	0.55	0.560	0.56	0.532
FCM-3-5	0.52	0.52	0.56	0.539	0.534	0.579	0.573		0.55	0.552	0.56	0.541
FCM-3-6	0.53	0.51	0.58	0.527	0.529	0.558	0.571	0.555	0.55	0.556	0.55	0.524
FCM-3-7	0.53	0.51	0.56	0.533	0.528	0.544	0.579	0.556	0.55	0.553	0.55	0.535
FCM-3-8	0.54	0.52	0.55	0.539	0.535	0.540	0.576	0.555	0.54	0.564	0.54	0.529
FCM-3-9	0.52	0.51	0.56	0.547	0.531	0.554	0.570	0.554	0.54	0.558	0.54	0.525
FCM-3-10	0.53	0.51	0.55	0.549	0.530	0.555	0.577	0.554	0.54	0.556	0.55	0.526
Mean	0.528	0.512	0.560	0.543	0.531	0.553	0.570	0.556	0.546	0.556	0.547	0.529
Std. Devn.	0.0063	0.0063	0.0094	0.0094	0.0024	0.0116	0.0076	0.0022	0.0052	0.0054	0.0095	0.0056
% RSD	1.20	1.22	1.68	1.74	0.45	2.10	1.34	0.39	0.95	0.98	1.73	1.06

STANDARD REFERENCE MATERIAL CDN-FCM-3

Participating Laboratories:

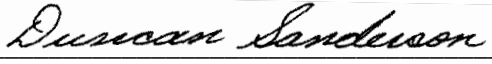
(not in same order as listed in table of results)

Acme Analytical Laboratories Ltd., Vancouver
Assayers Canada Ltd., Vancouver
Alex Stewart (Assayers) Argentina
ALS Chemex Laboratories, North Vancouver
Eco-Tech Laboratories Ltd., Kamloops
Genalysis Laboratory Services Ltd., Perth
GTK Laboratory, Finland
International Plasma Laboratories Ltd., Vancouver
OMAC Laboratory Ltd., Ireland
SGS-XRAL, Toronto
Teck Cominco - Global Discovery Laboratory, Vancouver
TSL Laboratories Ltd., Saskatoon


Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. or Barry Smee accept no liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by


Duncan Sanderson, Certified Assayer of B.C.

Geochemist


Dr. Barry Smee, Ph.D., P. Geo.