

CDN Resource Laboratories Ltd.

#2, 20148 – 102nd Avenue, Langley, B.C., Canada, V1M 4B4, 604-882-8422, Fax: 604-882-8466 (www.cdnlabs.com)

REFERENCE MATERIAL: CDN-GS-4C

Recommended value and the "Between Laboratory" two standard deviations

Gold concentration: 4.26 ± 0.22 g/t (30g Fire Assay / ICP)

Gold concentration: 4.25 ± 0.20 g/t (30g Fire Assay / Grav.)

PREPARED BY: CDN Resource Laboratories Ltd.
CERTIFIED BY: Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia
INDEPENDENT GEOCHEMIST: Dr. Barry Smee., Ph.D., P. Geo.
DATE OF CERTIFICATION: March 10, 2010

ORIGIN OF REFERENCE MATERIAL:

Standard CDN-GS-4C was prepared using ore supplied by Barrick Gold Inc. from their Goldstrike Mine in Nevada, USA. It is Carlin Style Mineralization in the prolific Northern Carlin Trend in Northern Nevada, USA. The source material is from Devonian carbonates of the Popovich Formation. Gold is strongly associated with pyrite and other sulfides including the arsenic minerals orpiment and realgar.

METHOD OF PREPARATION:

Reject ore material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 material was mixed for 5 days in a double-cone blender. Splits were taken and sent to 15 commercial laboratories for round robin assaying. Round robin results are displayed below:

ICP	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
CDN-GS-4C-1	4.23	4.18	4.52	4.16	4.11	4.01	4.46	4.20	4.34	4.41	4.34	4.10	4.26	4.23	4.21
CDN-GS-4C-2	4.48	4.20	4.44	4.13	4.13	4.08	4.56	4.30	4.29	4.40	4.21	4.19	4.22	4.47	4.15
CDN-GS-4C-3	4.23	4.16	4.37	4.20	4.21	4.05	4.47	4.27	4.32	4.39	4.30	4.18	4.29	4.43	4.09
CDN-GS-4C-4	4.24	4.20	4.27	4.25	4.22	4.13	4.58	4.20	4.38	4.44	4.30	4.16	4.24	4.41	4.15
CDN-GS-4C-5	4.46	4.27	4.59	4.14	4.21	4.09	4.47	4.29	4.28	4.38	4.25	4.11	4.27	4.39	4.17
CDN-GS-4C-6	4.18	4.34	4.60	4.26	4.27	4.07	4.35	4.36	4.26	4.33	4.32	4.13	4.22	4.39	4.11
CDN-GS-4C-7	4.27	4.26	4.46	4.30	4.24	4.06	4.40	4.28	4.29	4.34	4.27	4.19	4.30	4.25	4.20
CDN-GS-4C-8	4.15	4.41	4.38	4.16	4.24	4.02	4.27	4.22	4.31	4.34	4.25	4.14	4.26	4.43	4.16
CDN-GS-4C-9	4.29	4.29	4.59	4.27	4.20	4.06	4.28	4.13	4.28	4.39	4.37	4.00	4.32	4.35	4.12
CDN-GS-4C-10	4.21	4.33	4.46	4.15	4.27	4.13	4.28	4.18	4.22	4.41	4.27	3.95	4.31	4.43	4.13
Mean	4.27	4.26	4.47	4.20	4.21	4.07	4.41	4.24	4.30	4.38	4.29	4.12	4.27	4.38	4.15
Std. Dev'n	0.1109	0.0804	0.1092	0.0625	0.0533	0.0420	0.1146	0.0685	0.0440	0.0359	0.0476	0.0810	0.0349	0.0796	0.0381
%RSD	2.59	1.89	2.44	1.49	1.27	1.03	2.60	1.61	1.02	0.82	1.11	1.97	0.82	1.82	0.92
Gravimetric	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
CDN-GS-4C-1	4.32	4.16	4.80	4.23	4.33	4.146	4.26	4.31	4.42		4.3		4.19	4.120	4.30
CDN-GS-4C-2	4.32	4.23	4.53	4.26	4.33	4.065	4.23	4.24	4.39		4.2		4.23	4.310	4.33
CDN-GS-4C-3	4.35	4.15	4.33	4.30	4.36	4.058	4.20	4.25	4.51		4.1		4.23	4.467	4.27
CDN-GS-4C-4	4.28	4.26	4.57	4.33	4.34	4.121	4.13	4.22	4.56		4.2		4.22	4.067	4.30
CDN-GS-4C-5	4.26	4.25	4.77	4.26	4.30	4.012	4.24	4.30	4.22		4.1		4.24	4.061	4.33
CDN-GS-4C-6	4.24	4.11	4.33	4.20	4.32	4.091	4.26	4.33	4.37		4.4		4.24	4.223	4.20
CDN-GS-4C-7	4.23	4.19	4.50	4.33	4.26	4.097	4.15	4.23	4.58		4.2		4.19	4.061	4.30
CDN-GS-4C-8	4.27	4.26	4.67	4.31	4.31	4.146	4.28	4.24	4.38		4.3		4.20	4.122	4.27
CDN-GS-4C-9	4.08	4.18	4.30	4.23	4.26	4.103	4.28	4.23	4.43		4.5		4.22	4.238	4.23
CDN-GS-4C-10	4.25	4.26	4.67	4.19	4.33	4.198	4.26	4.17	4.42		4.4		4.20	4.000	4.23
Mean	4.26	4.21	4.55	4.26	4.31	4.10	4.23	4.25	4.43		4.27		4.22	4.17	4.28
Std. Dev'n	0.0742	0.0544	0.1832	0.0517	0.0327	0.0527	0.0528	0.0480	0.1042		0.1337		0.0196	0.1426	0.0443
%RSD	1.74	1.29	4.03	1.21	0.76	1.28	1.25	1.13	2.35		3.13		0.46	3.42	1.04

Note: Labs 10 and 12 did not report gravimetric results.
 Gravimetric data from Lab 3 was excluded for failing the "t" test.

REFERENCE MATERIAL: CDN-GS-4C

APPROXIMATE CHEMICAL COMPOSITION:

	Percent			Percent			ppm
SiO ₂	45.1		Na ₂ O	0.1		As	500
Al ₂ O ₃	5.2		MgO	8.9		Sb	60
Fe ₂ O ₃	2.7		K ₂ O	1.1			
CaO	14.3		TiO ₂	0.2			
MnO	0.1		LOI	20.8			
Total S	1.4		Total C	5.9			
Sulphide S	1.2		Inorganic C	5.3			

Statistical Procedures:

The final limits were calculated after first determining if all data was compatible within a spread normally expected for similar analytical methods done by reputable laboratories. Data from any one laboratory was removed from further calculations when the mean of all analyses from that laboratory failed a t test of the global means of the other laboratories. The mean and standard deviation were calculated using all remaining data. Any analysis that fell outside of the mean ± 2 standard deviations was removed from the ensuing data base. The mean and standard deviations were again calculated using the remaining data. This method is different from that used by Government agencies in that the actual "between-laboratory" standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.


Participating Laboratories: (not in same order as table of assays)

Acme Analytical Laboratories Ltd., Vancouver, B.C., Canada
Acme Analytical Laboratories Ltd., Santiago, Chile
Activation Laboratories, Ancaster, Ontario, Canada
Activation Laboratories, Thunder Bay, Ontario, Canada
Alaska Assay Lab, Anchorage, Alaska, USA
Alex Stewart Argentina SA
ALS Chemex, North Vancouver, B.C., Canada
ALS Chemex, Nevada, USA
Assayers Canada Ltd., Vancouver, B.C., Canada
Bourlamaque Assay Laboratories, Quebec, Canada
Eco Tech Laboratory Ltd., Kamloops, B.C., Canada
Labtium Inc., Finland
SGS, Toronto, Canada
SGS, Lima, Peru
TSL Laboratories Ltd., Saskatoon, SK, Canada

Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. nor Barry Smee accept any liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by

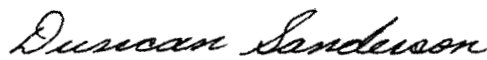

Duncan Sanderson, Certified Assayer of B.C.

Geochemist



Dr. Barry Smee, Ph.D., P. Geo.