

# CDN Resource Laboratories Ltd.

#2, 20148 – 102<sup>nd</sup> Ave, Langley, B.C., Canada, V1M 4B4, 604-882-8422, Fax: 604-882-8466 (www.cdnlabs.com)

## REFERENCE MATERIAL: CDN-ME-13

Recommended values and the “Between Lab” Two Standard Deviations

<i>Gold</i>	<i>0.148 g/t ± 0.024 g/t</i>	<i>provisional value (RSD = 10.4%)</i>
<i>Silver</i>	<i>76.5 g/t ± 6.8 g/t</i>	<i>Certified value</i>
<i>Copper</i>	<i>2.69 % ± 0.20 %</i>	<i>Certified value</i>
<i>Lead</i>	<i>1.70 % ± 0.09 %</i>	<i>Certified value</i>
<i>Zinc</i>	<i>18.48 % ± 1.83 %</i>	<i>Certified value</i>

**Note:** Standards with an RSD of near or less than 5% are certified; RSD's of between 5% and 15% are Provisional; RSD's over 15% are Indicated. Provisional and Indicated values cannot be used to monitor accuracy with a high degree of certainty.

**PREPARED BY:** CDN Resource Laboratories Ltd.  
**CERTIFIED BY:** Duncan Sanderson, B.Sc., Licensed Assayer of British Columbia  
**INDEPENDENT GEOCHEMIST:** Dr. Barry Smee., Ph.D., P. Geo.  
**DATE OF CERTIFICATION:** February 6, 2011

### **METHOD OF PREPARATION:**

Reject ore material was dried, crushed, pulverized and then passed through a 270 mesh screen. The +270 material was discarded. The -270 material was mixed for 5 days in a double-cone mixer. Splits were taken and sent to 15 laboratories for round robin assaying.

### **ORIGIN OF REFERENCE MATERIAL:**

The ore is described as massive to semi-massive sulphides from the Izok Lake orebody, an archean aged VMS deposit in the Slave structural province of Canada. It consists of pyrite, pyrrhotite, chalcopyrite, sphalerite and minor galena. Gangue minerals include quartz, chlorite, feldspar, cordierite, biotite, magnetite, anthophyllite and grunerite.

**Approximate chemical composition (from whole rock analysis) is as follows:**

	Percent		Percent
SiO <sub>2</sub>	11.4	MgO	1.4
Al <sub>2</sub> O <sub>3</sub>	1.4	K <sub>2</sub> O	0.2
Fe <sub>2</sub> O <sub>3</sub>	44.7	TiO <sub>2</sub>	<0.1
CaO	0.8	LOI	21.1
Na <sub>2</sub> O	0.5	S	35.1
C	0.1		

### **Statistical Procedures:**

The final limits were calculated after first determining if all data was compatible within a spread normally expected for similar analytical methods done by reputable laboratories. Data from any one laboratory was removed from further calculations when the mean of all analyses from that laboratory failed a t test of the global means of the other laboratories. The means and standard deviations were calculated using all remaining data. Any analysis that fell outside of the mean ±2 standard deviations was removed from the ensuing data base. The mean and standard deviations were again calculated using the remaining data. This method is different from that used by Government agencies in that the actual “between-laboratory” standard deviation is used in the calculations. This produces upper and lower limits that reflect actual individual analyses rather than a grouped set of analyses. The limits can therefore be used to monitor accuracy from individual analyses, unlike the Confidence Limits published on other standards.

### **Assay Procedures:**

**Au:** Fire assay pre-concentration, AA or ICP finish (10 or 15g sub-sample).  
**Ag, Cu, Pb, Zn:** 4-acid digestion, AA or ICP finish.

## REFERENCE MATERIAL CDN-ME-13

### Results from round-robin assaying:

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t	Au g/t
ME-13-1	0.15	0.122	0.16	0.134	0.20	0.165	0.157	0.155	0.157	0.136	0.140	0.15	0.146	0.12	0.19
ME-13-2	0.14	0.137	0.14	0.127	0.19	0.150	0.162	0.148	0.144	0.149	0.160	0.13	0.146	0.15	0.21
ME-13-3	0.16	0.155	0.13	0.146	0.19	0.160	0.161	0.153	0.138	0.160	0.140	0.14	0.135	0.13	0.21
ME-13-4	0.14	0.157	0.15	0.160	0.19	0.165	0.156	0.164	0.142	0.139	0.140	0.15	0.123	0.13	0.19
ME-13-5	0.14	0.156	0.15	0.163	0.19	0.170	0.148	0.133	0.156	0.138	0.140	0.15	0.163	0.15	0.21
ME-13-6	0.13	0.144	0.14	0.133	0.18	0.155	0.147	0.139	0.151	0.139	0.160	0.16	0.135	0.13	0.20
ME-13-7	0.14	0.162	0.15	0.122	0.23	0.155	0.147	0.166	0.140	0.143	0.150	0.17	0.142	0.13	0.22
ME-13-8	0.15	0.147	0.16	0.192	0.24	0.165	0.165	0.162	0.154	0.128	0.140	0.14	0.143	0.13	0.18
ME-13-9	0.14	0.146	0.17	0.163	0.19	0.170	0.160	0.153	0.157	0.143	0.140	0.15	0.112	0.15	0.19
ME-13-10	0.14	0.169	0.16	0.146	0.19	0.160	0.150	0.154	0.139	0.131	0.140	0.14	0.157	0.13	0.18
Mean	0.143	0.150	0.151	0.149	0.199	0.162	0.155	0.153	0.148	0.141	0.145	0.147	0.140	0.135	0.198
Std. Devn.	0.0082	0.0134	0.0120	0.0213	0.0197	0.0067	0.0068	0.0105	0.0079	0.0091	0.0085	0.0111	0.0150	0.0120	0.0140
% RSD	5.76	8.99	7.93	14.31	9.90	4.14	4.38	6.89	5.37	6.46	5.86	7.59	10.70	8.93	7.06
	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t	Ag g/t
ME-13-1	82	82	70.9	71.8	76.4	76.1	74.8	77	78	66.9	76.5	76	82	78.1	71.5
ME-13-2	80	78	71.2	72.1	75.2	76.7	77.9	78	78	71.9	77.0	77	83	80.4	63.5
ME-13-3	82	75	72.7	71.1	75.4	77.4	76.8	76	76	68.5	78.5	77	83	78.7	71.6
ME-13-4	81	80	71.6	69.6	72.4	77.6	77.6	79	77	68.3	74.5	77	83	78.4	71.1
ME-13-5	80	76	72.4	70.0	74.7	76.9	75.0	78	78	68.9	77.0	79	82	78.3	67.9
ME-13-6	81	78	72.7	69.6	75.8	77.3	75.9	78	76	70.0	77.0	78	83	78.7	70.4
ME-13-7	81	78	71.6	68.5	74.8	78.4	77.7	77	77	71.2	77.5	77	83	79.3	66.9
ME-13-8	83	78	71.9	70.0	76.0	76.2	75.5	77	76	69.7	78.5	79	82	77.7	67.7
ME-13-9	81	78	71.8	70.6	74.5	76.1	75.2	75	75	70.8	79.0	77	85	78.0	60.3
ME-13-10	82	78	71.7	69.6	76.1	77.2	76.8	75	77	73.2	76.5	78	86	80.5	72.3
Mean	81.3	78.1	71.9	70.3	75.1	77.0	76.3	77.0	76.8	69.9	77.2	77.3	83.1	78.8	68.3
Std. Devn.	0.9487	1.9120	0.5986	1.1110	1.1538	0.7430	1.1858	1.3333	1.0328	1.8840	1.2953	0.9559	1.3726	0.9954	3.9285
% RSD	1.17	2.45	0.83	1.58	1.54	0.97	1.55	1.73	1.34	2.69	1.68	1.24	1.65	1.26	5.75

**Note:** Au results from Laboratories 5 and 15 were removed for failing the “t” test.  
Ag results from Laboratory 15 were removed for failing the “t” test.

## REFERENCE MATERIAL CDN-ME-13

### Results from round-robin assaying:

	Lab 1	Lab 2	Lab 3	Lab 4	Lab 5	Lab 6	Lab 7	Lab 8	Lab 9	Lab 10	Lab 11	Lab 12	Lab 13	Lab 14	Lab 15
	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu	% Cu
ME-13-1	2.90	2.83	2.63	2.43	2.75	2.54	2.72	2.82	2.51	2.47	2.66	2.75	2.77	2.71	2.58
ME-13-2	2.83	2.76	2.64	2.42	2.73	2.51	2.81	2.86	2.73	2.49	2.65	2.75	2.82	2.70	2.61
ME-13-3	2.86	2.57	2.59	2.42	2.70	2.52	2.76	2.78	2.72	2.46	2.72	2.74	2.81	2.70	2.63
ME-13-4	2.86	2.73	2.63	2.39	2.67	2.52	2.79	2.82	2.60	2.53	2.75	2.74	2.80	2.70	2.62
ME-13-5	2.91	2.63	2.61	2.40	2.73	2.54	2.72	2.83	2.62	2.50	2.64	2.75	2.77	2.71	2.64
ME-13-6	2.85	2.70	2.65	2.45	2.68	2.55	2.75	2.79	2.67	2.52	2.60	2.72	2.80	2.72	2.58
ME-13-7	2.84	2.72	2.59	2.46	2.67	2.52	2.79	2.78	2.62	2.54	2.55	2.75	2.77	2.70	2.62
ME-13-8	2.87	2.71	2.61	2.49	2.77	2.54	2.75	2.78	2.67	2.43	2.69	2.72	2.77	2.72	2.62
ME-13-9	2.87	2.76	2.64	2.50	2.64	2.52	2.78	2.74	2.53	2.43	2.60	2.71	2.82	2.71	2.62
ME-13-10	2.84	2.73	2.60	2.48	2.71	2.53	2.76	2.72	2.62	2.44	2.73	2.75	2.81	2.71	2.53
Mean	2.862	2.714	2.619	2.444	2.705	2.529	2.762	2.792	2.629	2.481	2.659	2.738	2.795	2.708	2.605
Std. Devn.	0.0270	0.0717	0.0218	0.0381	0.0404	0.0129	0.0321	0.0421	0.0722	0.0412	0.0644	0.0155	0.0204	0.0068	0.0327
% RSD	0.94	2.64	0.83	1.56	1.49	0.51	1.16	1.51	2.75	1.66	2.42	0.57	0.73	0.25	1.26
	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb	% Pb
ME-13-1	1.80	1.72	1.54	1.75	1.70	1.60	1.66	1.71	1.66	1.68	1.67	1.69	1.70	1.74	1.70
ME-13-2	1.78	1.68	1.58	1.66	1.67	1.58	1.70	1.71	1.68	1.71	1.72	1.72	1.71	1.74	1.71
ME-13-3	1.81	1.61	1.57	1.71	1.71	1.59	1.71	1.70	1.63	1.69	1.76	1.69	1.74	1.75	1.73
ME-13-4	1.80	1.70	1.51	1.61	1.74	1.60	1.74	1.75	1.65	1.73	1.62	1.68	1.71	1.73	1.73
ME-13-5	1.80	1.65	1.52	1.70	1.81	1.60	1.68	1.72	1.67	1.72	1.71	1.69	1.70	1.73	1.71
ME-13-6	1.77	1.67	1.55	1.65	1.64	1.58	1.68	1.70	1.66	1.73	1.67	1.67	1.73	1.76	1.72
ME-13-7	1.78	1.68	1.57	1.66	1.78	1.58	1.72	1.66	1.66	1.74	1.68	1.69	1.72	1.77	1.75
ME-13-8	1.77	1.66	1.52	1.66	1.78	1.58	1.70	1.67	1.63	1.68	1.68	1.69	1.69	1.74	1.73
ME-13-9	1.79	1.68	1.56	1.66	1.77	1.58	1.69	1.66	1.62	1.67	1.68	1.68	1.71	1.74	1.73
ME-13-10	1.79	1.70	1.59	1.66	1.72	1.58	1.71	1.63	1.66	1.67	1.67	1.71	1.72	1.76	1.73
Mean	1.789	1.674	1.551	1.672	1.731	1.587	1.698	1.691	1.652	1.702	1.686	1.691	1.713	1.745	1.724
Std. Devn.	0.0137	0.0320	0.0277	0.0385	0.0543	0.0095	0.0226	0.0354	0.0193	0.0270	0.0372	0.0145	0.0140	0.0125	0.0143
% RSD	0.77	1.91	1.78	2.30	3.14	0.60	1.33	2.09	1.17	1.59	2.21	0.86	0.82	0.71	0.83
	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn	% Zn
ME-13-1	19.82	18.50	19.7	17.42	19.53	17.5	18.45	18.3	18.3	16.7	19.60	18.40	18.69	19.97	17.38
ME-13-2	19.38	18.10	19.9	17.67	19.57	17.4	18.67	18.2	18.6	16.7	19.40	18.40	18.72	19.93	17.39
ME-13-3	19.48	17.15	19.1	17.72	18.56	16.9	18.58	18.1	18.1	16.7	19.50	18.40	18.86	19.96	17.58
ME-13-4	19.54	18.25	19.2	17.50	18.98	17.2	18.29	18.9	17.8	17.1	19.50	18.40	18.64	19.85	17.52
ME-13-5	19.64	17.45	20.1	17.44	19.94	16.8	18.25	18.4	18.0	16.9	19.60	18.50	18.65	19.99	17.44
ME-13-6	19.38	17.70	19.9	17.89	20.13	17.5	18.23	18.5	18.4	17.1	19.10	18.30	18.86	19.67	17.50
ME-13-7	19.22	17.85	19.3	17.92	19.62	16.9	18.63	18.3	18.3	17.3	19.20	18.30	18.45	19.96	17.68
ME-13-8	19.57	17.85	19.4	18.17	20.36	16.9	18.45	17.9	17.8	16.5	19.40	18.20	18.74	19.60	17.66
ME-13-9	19.57	18.10	19.4	18.35	19.81	17.3	18.31	18.1	18.3	16.6	19.00	18.10	18.67	19.70	17.64
ME-13-10	19.31	18.15	19.2	17.98	19.52	17.4	18.55	17.9	18.3	16.7	19.50	18.40	18.72	19.70	17.65
Mean	19.49	17.91	19.52	17.81	19.60	17.18	18.44	18.26	18.19	16.83	19.38	18.34	18.70	19.83	17.54
Std. Devn.	0.1752	0.3992	0.3521	0.3129	0.5277	0.2781	0.1634	0.2989	0.2601	0.2584	0.2098	0.1174	0.1151	0.1482	0.1143
% RSD	0.90	2.23	1.80	1.76	2.69	1.62	0.89	1.64	1.43	1.54	1.08	0.64	0.62	0.75	0.65

**Note:** Cu results from Laboratory 4 were removed for failing the “t” test.  
Pb results from Laboratory 3 were removed for failing the “t” test.

**REFERENCE MATERIAL CDN-ME-13**

**Participating Laboratories:**

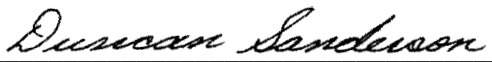
(not in same order as listed in table of results)

Acme Analytical Laboratories Ltd., Vancouver  
Actlabs-Ancaster, Ontario, Canada  
Actlabs-Thunder Bay, Ontario, Canada  
AGAT Laboratories, Ontario, Canada  
AHK Geochem, Alaska, USA  
ALS Chemex Laboratories, North Vancouver  
Genalysis Laboratory, Australia  
Labtium, Finland  
Omac Laboratories Ltd., Ireland  
Skyline Assayers and Laboratories, Arizona, USA  
SGS – Vancouver, B.C., Canada  
Stewart Group, Kamloops, B.C., Canada  
Alex Stewart Argentina SA  
TSL Laboratories Ltd., Saskatoon  
Ultra Trace Analytical Laboratories, Australia


Legal Notice:

This certificate and the reference material described in it have been prepared with due care and attention. However CDN Resource Laboratories Ltd. or Barry Smee accept no liability for any decisions or actions taken following the use of the reference material. Our liability is limited solely to the cost of the reference material.

Certified by

  
\_\_\_\_\_  
Duncan Sanderson, Certified Assayer of B.C.

Geochemist

  
\_\_\_\_\_  
Dr. Barry Smee, Ph.D., P. Geo.